

Door Connectors

Through the door – quick and easy



Primary application

Home care
Nursing homes
Hospital



Patient transfer from room to room can be a very difficult task since most overhead track systems require under door transfer. Ergolet has now introduced a unique door connector system for fast, easy, comfortable and, above all, safe transfers from room to room.

This door connector does not require costly changes to the electrical system since it is totally manual in operation. Simply pulling on the operating cord opens the connection from the one system to the other. To lock, pull on the cord again to open the locking mechanism, and click the traverse into place.

The locking mechanism allows both sets of track to connect smoothly and easily without the need to first line the traverse track up.

Safe mechanical systems ensure that the track cart stays in place in the track, even if the door connector fails. With its very compact size, Ergolet's door connector makes for a discreet, user friendly and functional overhead tracking unit.



Key standard functions

Compact and discreet
Flexible for fixed and travers systems
Manual operation


elite
healthcare

PRODUCT HIGHLIGHTS

- Manual operation – no costly electrical installations required
- Traverse clicks in place – no need for lining up
- Compact and discreet
- Several models to suit individual configurations
- Fail safe mechanical features




Traverse 135 mm
Connector 135 mm
Traverse 135 mm

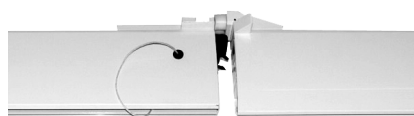

Fixed 135 mm
Traverse 135 mm




Traverse 135 mm
Connector 135 mm
Traverse 75 mm

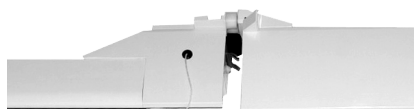

Fixed 75 mm
Traverse 135 mm



Traverse 75 mm
Connector 135 mm
Traverse 75 mm

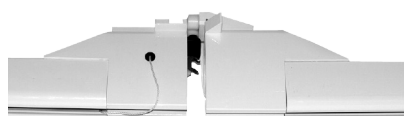


Fixed 75 mm
Traverse 75 mm


Version	Item-number
Traverse 135 mm - Connector 135 mm - Traverse 135 mm	016-01080
Fixed 135 mm - Traverse 135 mm	016-01081
Traverse 135 mm - Connector 135 mm - Traverse 75 mm	016-01072
Fixed 75 mm - Traverse 135 mm	016-01073
Traverse 75 mm - Connector 135 mm - Traverse 75 mm	016-01079
Fixed 75 mm - Traverse 75 mm	016-01078

TECHNICAL SPECIFICATION


	Connector 135mm	Traverse 75mm	Fixed 75mm
Length	300 mm	220 mm	220 mm
Height	145 mm	168 mm	145 mm
Width	100 mm	130 mm	100 mm

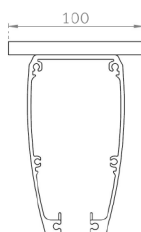




Drawings refer to Connector 135 mm

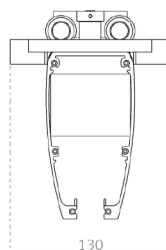

Drawings refer to Traverse 75 mm

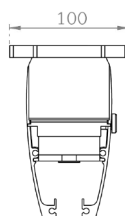

Drawings refer to Fixed 75 mm

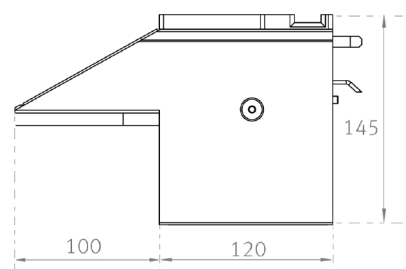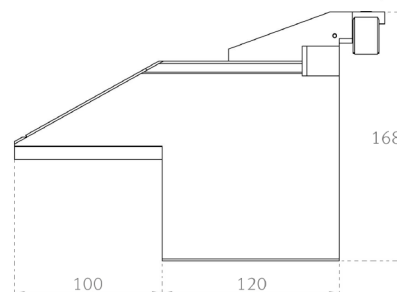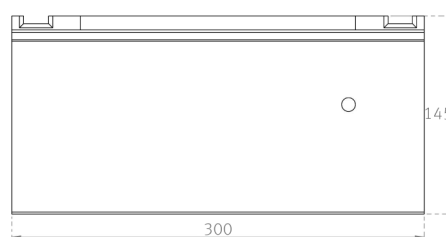

Consult the Technical manual for further specifications. All figures and measurements are for guidance only and may differ in actual situations. Ergolet reserves the right to changes without notice.

Product photos may include accessories. Consult your dealer.



date 06.01.2012 - ver. 3

Elite Healthcare Ltd.
Phone: +353 (0) 1 8903932
Fax: +353 (0) 1 8903937

Web: www.elitehealthcare.ie
Email: info@elitehealthcare.ie

Dublin: F1 Swords Enterprise Park,
Feltrim Road, Swords, Co Dublin,
K67 E984, Ireland.

Cork: 602 Harbour Point Business
Park, Little Island, Cork, T45 X624,
Ireland.

