



## UNIVERSAL 3-POINT with click-clips and head support

FOR THE POWERED POSITION SPREADER BAR (PPS) FOR LUNA® X5 CEILING LIFT

### USAGE

Used together with the PPS, this sling offers greater versatility for positioning and adjusting the patient during transfers. It allows movement from a semi-lying posture to an upright one in a smooth, controlled motion, and back again when needed.

### FEATURES

- The sling helps transferring the user from a wheelchair/chair or from lying position from a bed.
- The sling provides perfect support around the thighs and the high back, supporting the upper body and giving the user a feeling of comfort and safety.
- The clip is easy to attach by simply pulling it down over the spreader bar's stud.
- The clip solution is also easy to remove and gentle on fingers and caregivers, as it features a small release mechanism that is activated with a light press on the side of the clip, making detachment simple and comfortable.
- This locking mechanism also serves as an additional safety feature that prevents wear on the clip during repeated attachment and removal, reducing the risk of the clip becoming loose and unintentionally detaching.
- Location pocket for better sling fitting when sitting in a wheelchair.
- Straps on both sides of the sling to help the carer moving the user.
- Used with this sling, the PPS helps place the user correctly in the chair during transfers from bed to chair, including adjustments such as creating a more open seat angle.

### SPECIFICATIONS

SWL: 275 kg

**Material:**  
Polyester

**Sizes:** S-M-L-XL

**ITEM NUMBERS:**  
Universal Basic Clip  
with head support, polyester:

	Item no.
<b>S</b>	009-10504
<b>M</b>	009-10505
<b>L</b>	009-10506
<b>XL</b>	009-10507



**NEW!** We have introduced a new feature by equipping each strap of the clip sling with two clips, providing excellent adjustment possibilities.




## SLING SIZE GUIDE

The sling has a coloured strap on upper back with size and midline mark, which ensures correct placement of the sling. Colour code is presented on the size chart.

The image below shows a yellow colour code, which means that the sling is in size M (Medium):

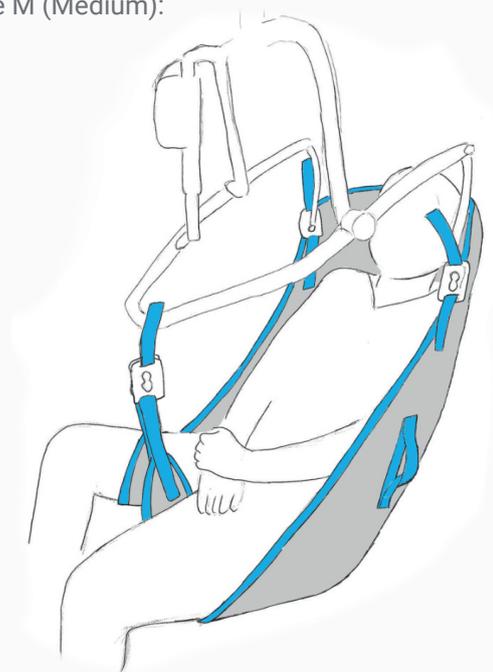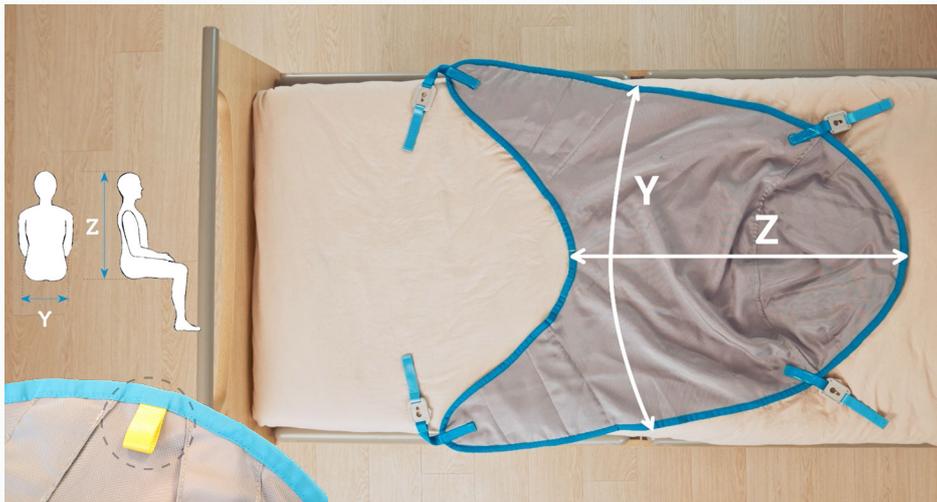

	S	M	L	XL
<b>Y</b> (hip measurement)	42 cm	48 cm	52 cm	68 cm
<b>Z</b> (measurement from head to lower back)	80 cm	91 cm	100 cm	100 cm

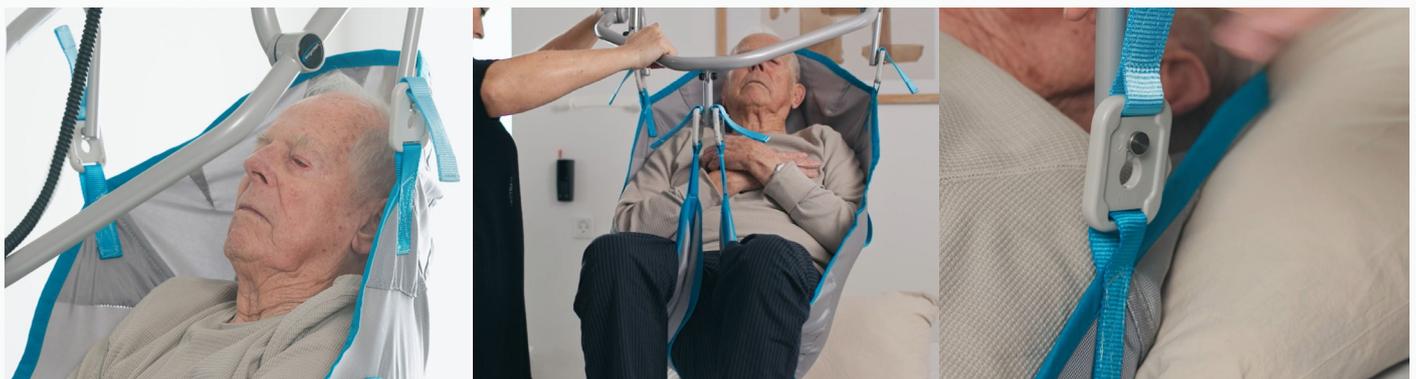

USE PROTECTIVE WASHING BAG

Please consult the technical manual for further specifications. All figures and measurements are for guidance only and may differ in actual situations. Winnicare Nordic reserves the right to changes without notice. Product photos may include accessories. Please consult your dealer for further information.

**elite**  
healthcare

20 Primeside Park  
Northwest Business Park  
Ballycoolin, Dublin 15. D15 E393

T: 01 8903932  
E: info@elitehealthcare.ie  
www.elitehealthcare.ie

